# C CARE

## KTS 210

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---

Service adres:  
Cresta International B.V  
Rolbrugweg 4 | NL-1332 AS ALMERE  
info@cresta-international.com  
www.cresta-international.com

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## INSTRUCTIONS

We congratulate you on purchasing our C-CARE electric heating pad KTS210 and thank you for the confidence you have placed in our brand and products.

We are convinced that with this purchase you have acquired a high quality product that you will enjoy using for many years to come.

Please read these operating instructions carefully before using the product, even if you are familiar with the type of product. It gives you useful tips to get the most out of your product.

If, after reading this manual, you still have questions about this product and its use, please send an e-mail to the address below; we will answer your question immediately.

CRESTA INTERNATIONAL B.V.

ROLBRUGWEG 4

NL-1332 AS ALMERE

[info@cresta-international.com](mailto:info@cresta-international.com)

[www.cresta-care.com](http://www.cresta-care.com)

Tel : + 31 36 3030063

# C CARE

## DESCRIPTION OF THE ICONS:

NOT TO BE USED WITH VERY YOUNG CHILDREN  
(0-3 YEARS)

READ THE INSTRUCTIONS

VOLTAGE


WATTAGE

OBSERVE ALL WARNINGS

FREQUENCY

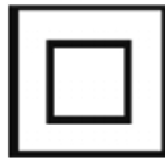PROTECTIONCLASS II

HOT SURFACE

KEEP CHILDREN AWAY FROM THE  
PRODUCT

DANGER FOR ELECTRIC  
SHOCK IF CABLE OR PLUG  
DAMAGED

ONLY USE INDOOR

# C CARE

DO NOT INSERT PINS

DO NOT USE FOLDED OR RUCKED

ATTENTION: ELECTRIC SHOCK

COVER MACHINE WASHABLE  
MAXIMUM WASH TEMPERATURE 30°C  
MILD PROCESS

DO NOT TUMBLE DRY

DO NOT BLEACH

DO NOT IRON

DO NOT DRY CLEAN

READ THIS MANUAL CAREFULLY AND FOLLOW THE INSTRUCTIONS

- This product is for indoor household use only, not for professional use in hospitals
- Only use this heating pad with the controller SS19
- The device is only for private use and not for use in hospitals
- Do not use the device folded or creased
- Only use the heating pad to warm your body. It gives you the heat in the places you want on your body.

# C CARE

- Only connect it to a network with the voltage stated on the care label
- The electric heating pad must not be used folded or wrinkled
- Do not use the heating pad on persons who cannot move without assistance, children under 3 years who cannot react to overheating, persons who are insensitive to heat or animals
- Persons who are only sensitive to heat up to a certain extend may only use the heating pad under supervision
- The appliance is not to be used by children over the age of three unless the controls have been pre-set by a parent or guardian, and unless the child has been adequately instructed on how to operate the controls safely.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory and mental capabilities or lack of experience and knowledge can use this product if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved.
- Children should not play with the product
- Cleaning and user maintenance shall not be made by children without supervision.
- The product may be used by children between 3 and 8 years of age under supervision and with the controller always set to the lowest temperature setting
- The cables and controller of the device can present a danger of entanglement, strangulation, tripping or stepping on if placed incorrectly.  
The user must ensure that excess cable ties and electrical cables are securely installed.
- Do not pierce with, not attach safety pins or other spiky or sharp objects to the heating pad.
- Do not pull or twist the power cord
- Do not switch the product by placing a timer between it
- Do not expose the switch to moisture or liquids
- If the fuse is blown or if the thermal connection is broken, return the heating pad to the store where you bought it.
- If the power cord is damaged, it must be replaced by the importer or its service department, or by a qualified person in order to avoid a hazardous situation.

# C CARE

- Do not place heavy objects such as a suitcase on the heating pad when it is switched on.
- Do not place any heat sources such as a hot water bottle on the heating pad when it is switched on
- Check the heating pad regularly for defects or signs of wear. If such signs appear or the heating pad has not been used properly, discontinue use and return it to the store where you purchased it.
- A damaged cord or plug can cause an electric shock
- Make sure that you do not fall asleep with the heating pad. Prolonged use can lead to skin burns
- If the device is used for an extended period of time, for example when the user falls asleep, the controls should be set to position 1.
- Set the device to temperature setting 1 before falling asleep or using it for a long time.
- Do not use the heating pad on parts of your body that are burned, injured or swollen
- Do not use the heating pad when it is damp or wet and only use it in dry environments
- Make sure the switch is not covered. The switch gets hot during use and must be able to dissipate this heat
- Do not fold the pad when it is still warm
- The appliance is not to be used by young children over the age of three unless the controls have been pre-set by a parent or guardian, and unless the child has been adequately instructed on how to operate the controls safely.
- The manufacturer is not responsible for damage resulting from incorrect use.
- Make sure the cable does not run over sharp corners and make sure that the cable does not run over any other heat source

CHILDREN UNDER THE AGE OF THREE ARE NOT TO USE THIS APPLIANCE DUE TO THEIR INABILITY TO REACT TO OVERHEATING.

THIS APPLIANCE MUST NOT TO BE USED BY PERSONS INSENSITIVE TO HEAT AND OTHER VERY VULNERABLE PERSONS WHO ARE UNABLE TO REACT TO OVERHEATING

THE CABLES AND CONTROL UNIT OF THIS APPLIANCE CAN GIVE RISE TO ENTANGLEMENT, STRANGULATION, TRIPPING OR STEPPING ON IF PLACED INCORRECTLY. THE USER MUST ENSURE THAT EXCESS CABLES TIES AND ELECTRICAL CABLES ARE SECURELY INSTALLED

## DESCRIPTION OF THE PRODUCT

This heat pad consists of three components:

1. The heat pad
2. The protective cover
3. The temperature control unit

When using, preferably place the heating pad in the protective cover. The protective cover is machine washable at 30 ° C.

The cable with the temperature control unit can be detached from the heating pad for easier storage.

### ATTENTION:

The heating pad may give off a plastic odor when used for the first time. This scent is harmless and will disappear after some use.

## SWITCHING ON AND ADJUST TEMPERATURE

To use the product, first attach the cable with plug to the heating pad and then plug the plug into the wall socket.

### SWITCH ON:

Press the on / off button  on the controller. The display lights up with the number 1 and the number 90. Here, the number 1 stands for the temperature setting and the number 90 for the set timer. By pressing the + or - you adjust the temperature setting to higher or lower. You can select the temperature settings from 1 to 6. For long-term use, we recommend using a maximum of setting 3

### SWITCH OFF:

Press the on / off button  on the controller to turn off the heating pad

## TEMPERATURE SELECTION

Select the temperature setting that is comfortable for you. For rapid heating, we recommend briefly using setting 6

## TEMPERATURE SETTINGS

Position 0	:	The product is switched off
Settings 1 and 2	:	Low heat
Settings 3 and 4	:	Average heat
Settings 5 and 6	:	Maximum heat

## AUTOMATIC POWER OFF

The heating pad is equipped with an adjustable timer that automatically switches the product off after 30, 60 or 90 minutes. When turned on, the default timer setting is 90 minutes.

Press the clock button  on the controller to change the timer setting to 30 or 60 minutes. If you still want to use the product after the automatic switch-off, switch the cushion on again as mentioned above. Please note that the temperature setting is 1 again and the timer is set to the standard setting of 90 minutes.

If you no longer want to use the heating pad, turn off the pad and remove the plug from the socket. If you want to store the heating pad for a long time, first remove the plug from the socket and then disconnect the cable with switch from the heating pad for easier storage

## MAINTENANCE AND CLEANING

Use the heating pad together with the textile cover to prevent grease or damage from using lotions or oils.

**CAUTION:** To avoid electric shock, before cleaning the product, remove the plug from the socket. Then let the pillow cool for at least 10 minutes.

Also disconnect the cable with controller from the heating pad

You can wash the heating pad by hand and the cover in the washing machine on a delicate cycle of 30°C. Use a delicate detergent in the correct dosage as indicated by the manufacturer. Small stains can be removed with a damp cloth or sponge and a mild soap solution. Do not use bleach for cleaning. Rinse several times to thoroughly remove all residual detergent. Do not use the heating pad again until it is completely dry.

Do not use chemicals or abrasives to clean the product, as this may damage the product.

The product cannot be chemically cleaned. The protective cover can be dry cleaned or cleaned according to the instructions on the product label

Never immerse the switch in water or other liquid. This can damage the product.

Do not hang the pillow to dry with clothespins, but lay it flat on a moisture-absorbing surface such as a large towel.

# C CARE

The product cannot be tumble dried. Also, do not turn on the product to dry it or iron it. Do not use the product again until it is completely dry.


## STORAGE

Let the heating pad cool down before putting it away. To avoid damage, do not place heavy objects on the heating pad.

If you do not want to use the heating pad for a longer period of time, place it back in the original packaging so that it is optimally protected.

## ENVIRONMENT FRIENDLY DISPOSAL

The European Directive 2002/96/EC on Waste Electrical and Electronic Equipment – WEEE – requires that old household electrical appliances should not be disposed of in the normal unsorted municipal waste stream. Old appliances should be collected separately in order to optimize the recovery and recycling of the materials they contain and reduce the impact on human health and the environment.

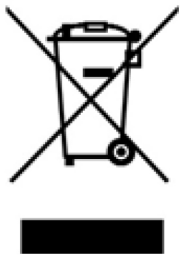The crossed out “wheeled bin” symbol on the product reminds you of your obligation that when you dispose of the appliance, it should be separately collected. Contact your local authority for information concerning the correct disposal of your old appliance.

## WARRANTY

CRESTA INTERNATIONAL B.V. warrants this product to be free from manufacturing 'defects and material errors for a period of 36 months from the date of purchase. In the event of a defect within this period, CRESTA will at its discretion repair or replace the device free of charge for the same or an equivalent device. This warranty does not apply in the event of normal wear and tear or if, in the opinion of CRESTA INTERNATIONAL B.V. the defect was caused by improper use, falling or impact, water damage, leaking batteries, if applicable, and the like. The warranty is only applicable if the warranty claim can be demonstrated with the purchase receipt.

---

Do not try to open the product and repair it yourself. Your warranty claims will lapse

---

### Requirements:

If a repair is necessary, you can return the product to the store where you bought it. You can also send the product directly to us. Send a message in advance with your complaint description and do not forget to state your full address details - not a PO box - when sending so that we can return the product to you after repair or replacement. If you send the product to us, you must frank it sufficiently. After repair or exchange you will receive your product from us delivered free.

### Service address:

Cresta International B.V

Rolbrugweg 4 | NL-1332 AS ALMERE

info@cresta-international.com

[www.cresta-international.com](http://www.cresta-international.com)

Printing errors and modifications reserved.

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