
C CARE

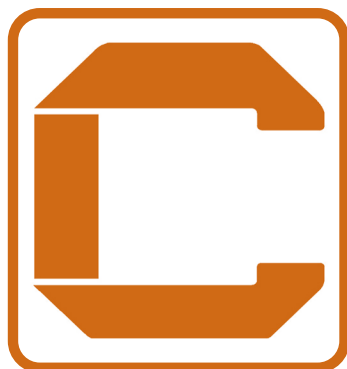KTS 220

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Service adres:

Cresta International B.V
Rolbrugweg 4 | NL-1332 AS ALMERE
info@cresta-international.com
www.cresta-international.com

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INSTRUCTIONS

Congratulations on your purchase of the C Care KTS220 Neck & Shoulder Heating Pad and thank you for the confidence you have placed in our brand and products.

We are convinced that with this purchase you have bought a high-quality product that you will enjoy for years to come.

Read this manual carefully before using the product, even if you are familiar with the type of product. It gives you useful tips to get the most out of your product.

If, after reading this manual, you have any questions about this product and its use, please send an e-mail to the address below; we will answer your question immediately:

CRESTA INTERNATIONAL BV

Rolbrugweg 4

NL-1332 AS ALMERE

Nederland

info@cresta-international.com

www.cresta-care.com

Tel: +31 36 3030063

C CARE

IMPORTANT INFORMATION, KEEP THIS FOR POSSIBLE FUTURE REFERENCE,

Please read the instruction manual carefully before using this product. Pay particular attention to the safety regulations. Keep this manual in an easily accessible place. If you hand over the product to someone, please also give this manual.

Read the instructions

Do not use if folded or wrinkled

Do not insert pins

Do not use with young children
(0-3 year)

READ THIS MANUAL CAREFULLY AND FOLLOW THE INSTRUCTIONS

- This product is for indoor household use only, not for professional use in hospitals
- Only use this warming jacket with the controller SS19
- Only use the warming jacket to warm your body. It gives you the heat in the places you want on your body.
- Only connect it to a network with the voltage stated on the care label
- The electric warming jacket must not be used folded or wrinkled

C CARE

- Do not use the warming jacket on persons who cannot move without assistance, children under 3 years who cannot react to overheating, persons who are insensitive to heat or animals
- Persons who are only sensitive to heat up to a certain extend may only use the heating pad under supervision
- The appliance is not to be used by children over the age of three unless the controls have been pre-set by a parent or guardian, and unless the child has been adequately instructed on how to operate the controls safely.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory and mental capabilities or lack of experience and knowledge can use this product if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved.
- Children should not play with the product
- Cleaning and user maintenance shall not be made by children without supervision.
- The product may be used by children between 3 and 8 years of age under supervision and with the controller always set to the lowest temperature setting
- The cables and controller of the device can present a danger of entanglement, strangulation, tripping or stepping on if placed incorrectly.
The user must ensure that excess cable ties and electrical cables are securely installed.
- Do not pierce with, nor attach safety pins or other spiky or sharp objects to the warming jacket.
- Do not pull or twist the power cord
- Do not switch the product by placing a timer between the product and the socket
- Do not expose the switch to moisture or liquids
- If the fuse is blown or if the thermal connection is broken, return the warming jacket to the store where you bought it.
- If the power cord is damaged, it must be replaced by the importer or its service department, or by a qualified person in order to avoid a hazardous situation.
- Do not place heavy objects such as a suitcase on the warming jacket when it is switched on.

C CARE

- Do not place any heat sources such as a hot water bottle on the warming jacket when it is switched on
- Check the warming jacket regularly for defects or signs of wear. If such signs appear or the heating pad has not been used properly, discontinue use and return it to the store where you purchased it.
- A damaged cord or plug can cause an electric shock
- Make sure that you do not fall asleep with the warming jacket. Prolonged use can lead to skin burns
- If the device is used for an extended period of time, for example when the user falls asleep, the controls should be set to position 1.
- Set the device to temperature setting 1 before falling asleep or using it for a long time.
- Do not use the warming jacket on parts of your body that are burned, injured or swollen
- Do not use the warming jacket when it is damp or wet and only use it in dry environments
- Make sure the switch is not covered. The switch gets hot during use and must be able to dissipate this heat
- Do not fold the jacket when it is still warm
- The appliance is not to be used by young children over the age of three unless the controls have been pre-set by a parent or guardian, and unless the child has been adequately instructed on how to operate the controls safely.
- The manufacturer is not responsible for damage resulting from incorrect use.
- Make sure the cable does not run over sharp corners and make sure that the cable does not run over any other heat source

CHILDREN UNDER THE AGE OF THREE ARE NOT TO USE THIS APPLIANCE DUE TO THEIR INABILITY TO REACT TO OVERHEATING

THIS APPLIANCE MAY NOT BE USED BY PERSONS INSENSITIVE TO HEAT AND OTHER VERY VULNERABLE PERSONS WHO ARE UNABLE TO REACT TO OVERHEATING

THE CABLES AND CONTROL UNIT OF THIS APPLIANCE CAN GIVE RISE TO ENTANGLEMENT, STRANGULATION, TRIPPING OR

C CARE

STEPPING ON IF PLACED INCORRECTLY. THE USE MUST ENSURE THAT EXCESS CABLES TIES AND ELETRICAL CABLES ARE SECURELY INSTALLED.

Please follow the washing instructions and possible after-treatment carefully. The instructions below explain exactly what you can and cannot do with the product

	DO NOT DRY CLEAN
	DO NOT BLEACH
	MAXIMUM WASHING TEMPERATURE 30°C - MILD PROCESS
	DO NOT TUMBLE DRY
	DO NOT IRON

PRODUCT APPLICATION

The KTS220 has been specially developed to provide heat locally around the neck, shoulders and back to loosen stiff muscles or relieve pain.

When storing the appliance allow it to cool down before folding

Do not crease the appliance by placing items on top of it before storage.

PRODUCT VIEW

Heating pad, cable and controller

USE OF THE PRODUCT

Apply the heating pad around the neck and close the snaps. Secure the cushion around your waist with the elastic cord with buckle. You can make the elastic wider or shorter for the correct fit.

Plug in the power cord and turn on the pad by pressing the on / off button  on the controller. The numbers 1 and 90 appear in the display. The number 1 stands for the temperature setting, the number 90 for the set timer.

You can adjust the temperature by pressing the + and – buttons on the controller

Position 0	:	product is switched off
Position 1 + 2	:	minimum heat
Position 3 + 4	:	medium heat
Position 5 + 6	:	maximum heat

In order to heat up the cushion quickly, we advise you to use position 6. If you want to use the product for a longer period of time, set position 3.

Remove the plug from the wall socket after use.

Attention:

This product is designed for home use, not for professional use.

Do not use this product on babies, infants, persons insensitive to heat, unsupervised disabled people and pets.

AUTOMATIC SWITCH OFF

The heating pad is equipped with an adjustable timer that automatically switches the product off after 30, 60 or 90 minutes. When switched on, the default setting for the timer is 90 minutes.

C CARE

Press the clock button  on the controller to change the timer setting to 30 or 60 minutes. If you want to continue using the product after it has turned off automatically, turn the pillow back on as described above. Please note that the temperature setting is back to 1 and the timer is set to the standard setting of 90 minutes.

When you no longer want to use the heating pad, switch the pad off and unplug it from the socket. If you want to store the heating pad for a longer period of time, first unplug the plug from the socket and then pull the switch off the cable from the heating pad for easier storage

CLEANING AND MAINTENANCE

Before cleaning the product, unplug the power cord and unplug the cable from the pad to prevent the controller from getting wet.

Remove small stains with a soft cloth or damp sponge with a few drops of mild detergent.

The cushion is not machine washable, cannot be tumble-dried, may not be wrung out, may not be dry cleaned and may not be ironed.

Do not use the pad until it is completely dry.

Limit the number of cleanings to avoid extra wear.

SPECIFICATIONS

Product	:	Heating pad for neck and shoulders
Model	:	KTS220
Consumption	:	100 W

ENVIRONMENT FRIENDLY DISPOSAL

The European Directive 2002/96/EC on Waste Electrical and Electronic Equipment - WEEE - requires that old household electrical appliances should not be disposed of in the normal unsorted municipal waste stream. Old appliances should be collected separately in order to optimize the recovery and recycling of the materials they contain and reduce the impact on human health and the environment.

The crossed out “wheeled bin” symbol on the product reminds you of your obligation that when you dispose of the appliance, it should be separately collected. Contact your local authority for information concerning the correct disposal of your old appliance.

WARRANTY

CRESTA INTERNATIONAL B.V. warrants this product to be free from manufacturing 'defects and material errors for a period of 36 months from the date of purchase. In the event of a defect within this period, CRESTA will at its discretion repair or replace the device free of charge for the same or an equivalent device. This warranty does not apply in the event of normal wear and tear or if, in the opinion of CRESTA INTERNATIONAL B.V. the defect was caused by improper use, falling or impact, water damage, leaking batteries, if applicable, and the like. The warranty is only applicable if the warranty claim can be demonstrated with the purchase receipt.

Do not try to open the product and repair it yourself. Your warranty claims will lapse

Requirements:

If a repair is necessary, you can return the product to the store where you bought it. You can also send the product directly to us. Send a message in advance with your complaint description and do not forget to state your full address details - not a PO box - when sending so that we can return the product to you after repair or replacement. If you send the product to us, you must frank it sufficiently. After repair or exchange you will receive your product from us delivered free.

Service address:

Cresta International B.V

Rolbrugweg 4 | NL-1332 AS ALMERE

info@cresta-international.com

www.cresta-international.com

Printing errors and modifications reserved

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